


CARDIOVIT AT-102 G2

Easy, smart and state-of-the-art

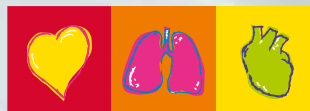

SCHILLER
The Art of Diagnostics

CARDIOVIT AT-102 G2

The new generation ECG: CARDIOVIT AT-102 G2

- ∴ Best signal quality and diagnostic algorithms for adult and paediatric patients
- ∴ Fast bidirectional Wi-Fi communication with HIS/EMR
- ∴ Spirometry and basic exercise ECG optionally available



SpiroScout SP plus, for spirometry tests with the CARDIOVIT AT-102 G2

STATE-OF-THE-ART TECHNOLOGY AND CONNECTIVITY

The CARDIOVIT AT-102 G2 features an 8" colour display with touch function keys. Its sealed alphanumeric keyboard ensures easy cleaning, making it highly suitable for the daily clinical routine. Smart battery technology guarantees more than 8 hours of ECG recording.

Bidirectional communication allows for easy data access and fast transmission of all reports to EMR/HIS systems in PDF, DICOM or HL7 format. Standard Wi-Fi and LAN connectivity, combined with ECG review on the display and a large memory, supports paperless workflows and cost saving. Customizable printouts are available on the internal high-resolution A4 thermal printer.

HIGH EFFICIENCY IN ANY WORKFLOW

ECG acquisition has never been easier and faster: Enter or download patient data with PDQ or worklist, acquire ECG, save, print and transmit data within twenty seconds using only one navigation key.




EASY 1-2-3 STEPS

- 1: Enter or download patient data
- 2: Check signal quality and perform ECG
- 3: Review ECG before storing, exporting or printing

THE CARDIOVIT AT-102 G2 DIAGNOSTIC TOOLS

- ∴ Resting ECG with the SCHILLER ECG analysis program ETM as well as ETM Sport, the automatic interpretation for athletes.
- ∴ Resting rhythm recording up to 10 minutes.
- ∴ Optimal signal quality for paediatric ECGs.
- ∴ Basic exercise ECG (optional) with a large selection of protocols for treadmills and ergometers.
- ∴ SCHILLER's Culprit Coronary Artery Algorithm (CCAA, optional) locates the obstruction site in the coronary artery and thereby the size of the cardiac area at risk.

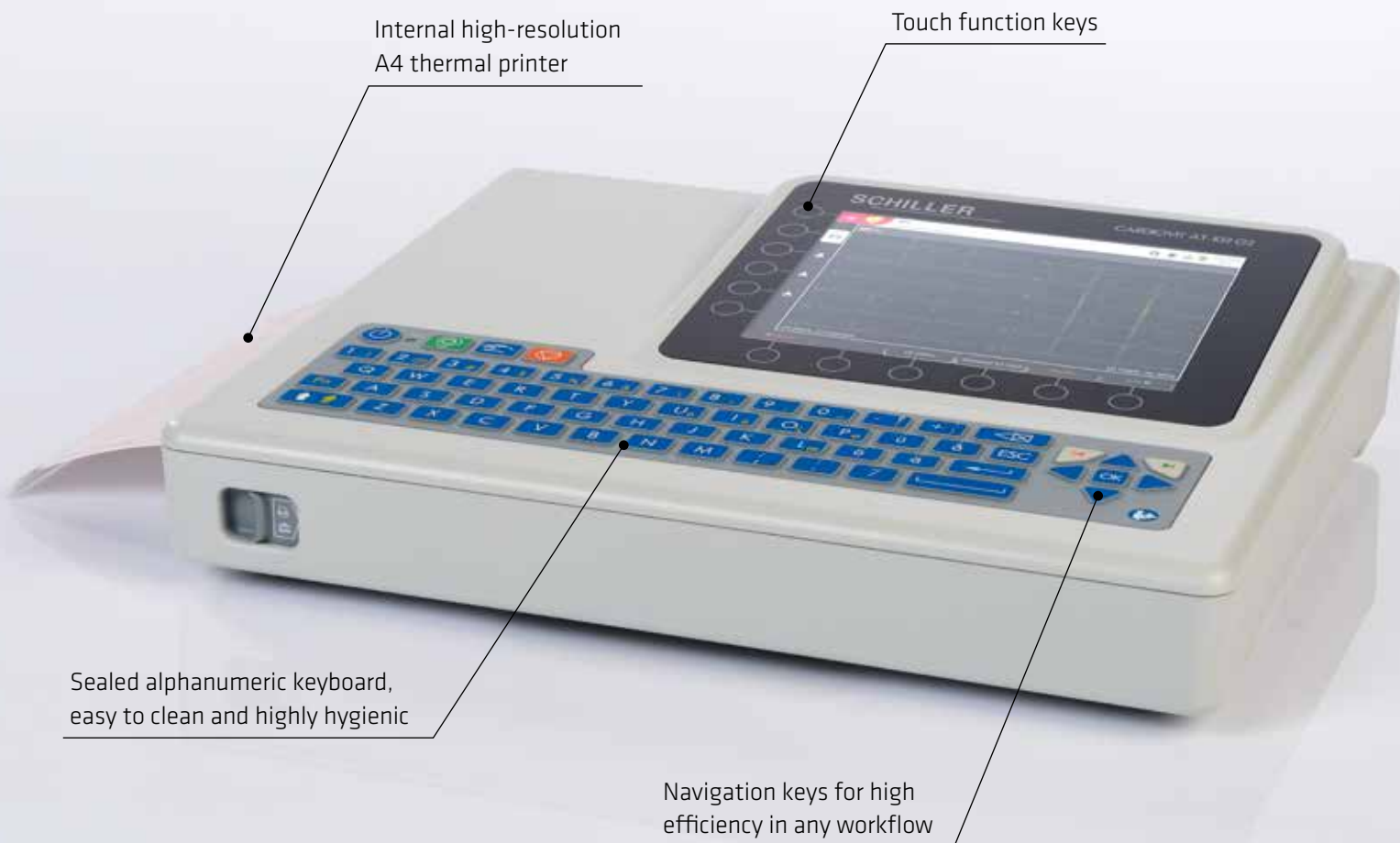

The CARDIOVIT AT-102 G2 is also available with a trolley.

ADD SPIROMETRY TO YOUR CARDIOVIT AT-102 G2

With the SpiroScout SP plus, the CARDIOVIT AT-102 G2 also becomes a spirometer. Thanks to GANSHORN's ultrasound technology, spirometry recordings are highly accurate and the sensor is maintenance- and calibration-free. The disposable ScoutTubes ensure maximum hygiene and patient safety. ScoutTubes are CO₂ neutral.

The spirometry option offers FVC, SVC and MVV with a large selection of reference standards, pre and post tests, as well as interpretation of the FVC results.

The device provides feedback on the procedure's quality during the test and after each trial.

Complete spirometry results can be printed on thermal paper as well as exported to EMR/HIS systems.



An incentive screen is included for trials with children.



Americas
SCHILLER Americas Inc.
Doral, Florida 33172

North America:
Phone +1 786 845 06 20
Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com

Latin America & Caribbean
Phone +1 305 591 11 21
Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com



Asia
SCHILLER Asia-Pacific / Malaysia
52200 Kuala Lumpur, Malaysia
Phone +603 6272 3033
sales@schiller.com.my
www.schiller-asia.com



Austria
SCHILLER Handelsgesellschaft m.b.H.
A-4040 Linz
Phone +43 732 70 99 0
Fax +43 732 757 000
sales@schiller.at
www.schiller.at



China
Alfred Schiller (Beijing) Medical Equipment
Co., Ltd.
100015 Beijing, China
Phone +86 010 52007020
Fax +86 010 52007020-8016
info@schillerchina.com
www.schillermedical.cn



Croatia
SCHILLER d.o.o.
10000 Zagreb
Phone +385 1 309 66 59
Fax +385 1 309 66 60
info@schillerzg.hr
www.schiller.hr



France
SCHILLER Médical
F-67160 Wissembourg
Phone +33 3 88 63 36 00
Fax +33 3 88 63 36 49
info@schiller.fr
www.schiller-medical.com



France (distribution France)
SCHILLER France S.A.S.
F-77608 Bussy St Georges
Phone +33 1 64 66 50 00
Fax +33 1 64 66 50 10
contact@schillerfrance.fr
www.schiller-france.com



Germany
SCHILLER Medizintechnik GmbH
D-85622 Feldkirchen b. München
Phone +49 89 62 99 81 0
Fax +49 89 609 50 90
info@schillermed.de
www.schillermed.de



India
SCHILLER Healthcare India Pvt. Ltd.
Mumbai - 400 059, India
Phone +91 22 6152 3333 / 2920 9141
Fax +91 22 2920 9142
sales@schillerindia.com
www.schillerindia.com



Poland
SCHILLER Poland Sp. z o.o.
PL-02-729 Warszawa
Phone +48 22 843 20 89 / +48 22 647 35 90
Fax +48 22 843 20 89
schiller@schiller.pl
www.schiller.pl



Russia & C.I.S.
AO SCHILLER.RU
125040 Moscow, Russia
Phone +7 (495) 970 11 33
Fax +7 (495) 956 29 10
mail@schiller.ru
www.schiller.ru



Serbia
SCHILLER d.o.o.
11010 Beograd
Phone +381 11 39 79 508
Fax +381 11 39 79 518
info@schiller.rs
www.schiller.rs



Slovenia
SCHILLER d.o.o.
2310 Slovenska Bistrica
Phone +386 2 843 00 56
Fax +386 2 843 00 57
info@schiller.si
www.schiller.si



Spain
SCHILLER ESPAÑA, S.A.
ES-28232 Las Rozas/Madrid
Phone +34 91 713 01 76
Fax +34 91 355 79 33
schiller@schiller.es
www.schiller.es



Switzerland
SCHILLER-Reomed AG
CH-8912 Obfelden
Phone +41 44 744 30 00
Fax +41 44 740 37 10
sales@schiller-reomed.ch
www.schiller-reomed.ch



Turkey
SCHILLER TÜRKIYE
Okmeydani-Sisli - Istanbul
Phone +90 212 210 8681 (pbx)
Fax +90 212 210 8684
info@schiller.com.tr
www.schiller-turkiye.com



United Kingdom
SCHILLER UK Ltd.
Coatbridge, ML5 2DA
Phone +44 1619 072 555
sales@schilleruk.com
www.schilleruk.com

Device availability in your market is subject to regulatory approval.

