
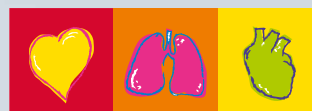


CARDIOVIT FT-1

Where performance meets portability


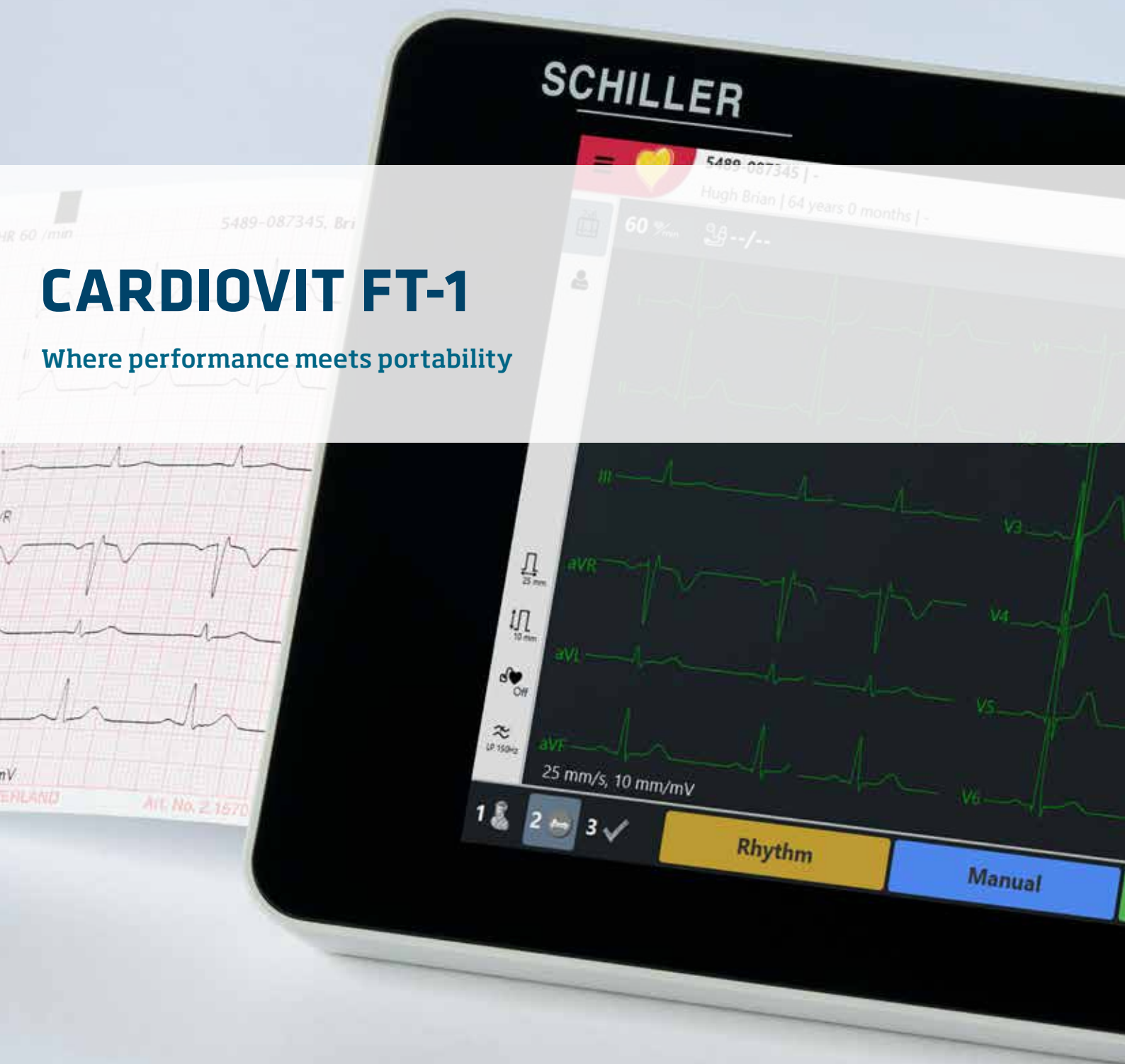


SCHILLER
The Art of Diagnostics

CARDIOVIT FT-1

The SCHILLER 12-lead ECG:

- ✧ Unique design
- ✧ Intuitive 3-step workflow for effortless operation
- ✧ Fast and secure bidirectional Wi-Fi communication
- ✧ Outstanding signal quality for adult and paediatric ECG

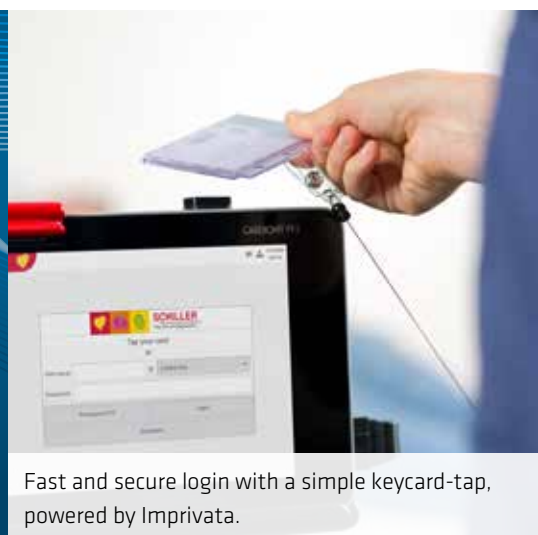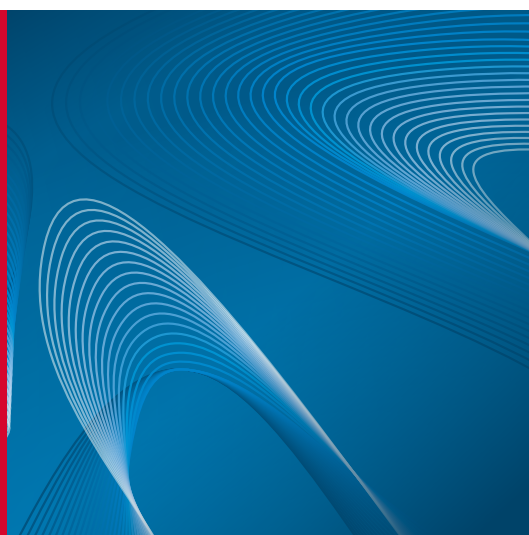

Fast and secure login with a simple keycard-tap, powered by Imprivata.

A UNIQUE DESIGN WITH INTEGRATED PRINTER

The CARDIOVIT FT-1 is an ultra-portable 12-lead ECG device, weighing just over 1 kg including the battery. It features an 8-inch high-resolution multi-touch screen and an integrated thermal printer.

EXTENSIVE CONNECTIVITY

Bidirectional communication enables easy data access via worklist or PDQ and fast transmission of all reports to EMR/HIS systems. Secure Wi-Fi connectivity, combined with on-screen ECG review and high storage capacity, enables reliable paperless workflows and reduces costs.

ADVANCED CYBERSECURITY

CARDIOVIT FT-1 offers exceptional cybersecurity through:

- ✧ Security-hardened Linux kernel to minimise the risk in case of cyberattacks
- ✧ Configurable access control with username, password, and privileges, as well as Imprivata's Enterprise Access Management (EAM) with tap-and-go access, to prevent unauthorised access
- ✧ Wi-Fi protocols including certificate-based authentication to ensure secure access to the hospital IT infrastructure
- ✧ Enhanced security by disabling unused USB ports

HIGH-END DIAGNOSTIC TOOLS

RESTING RHYTHM: CARDIOVIT FT-1 records all 12 leads for up to 4 minutes, reducing the risk of missing clinically relevant rhythm events. Recordings can be printed, stored, reviewed, and transmitted.

ETM WITH ETM SPORT (OPTIONAL): SCHILLER's renowned, simultaneous 12-lead ECG interpretation program for uncompromising quality and reliability. ETM Sport for the interpretation of athletes' hearts distinguishes between physiological and pathological ECG variations, thereby reducing the number of false-positive results.

CCAA (OPTIONAL): SCHILLER's Culprit Coronary Artery Algorithm is used in emergency departments and outside the hospital whenever an immediate decision on chest pain is required. It locates the obstruction site in the coronary artery and thereby determines the size of the cardiac area at risk.

INTEGRATED PRINTER


High-end integrated thermal printer

HOOK-UP ADVISER

The anatomical hook-up adviser helps to place the electrodes correctly, with zoom-in and rotate functions for an optimal view.

LEAD REVERSAL DETECTION

In case of lead reversal, misplaced electrodes are immediately highlighted before ECG acquisition.

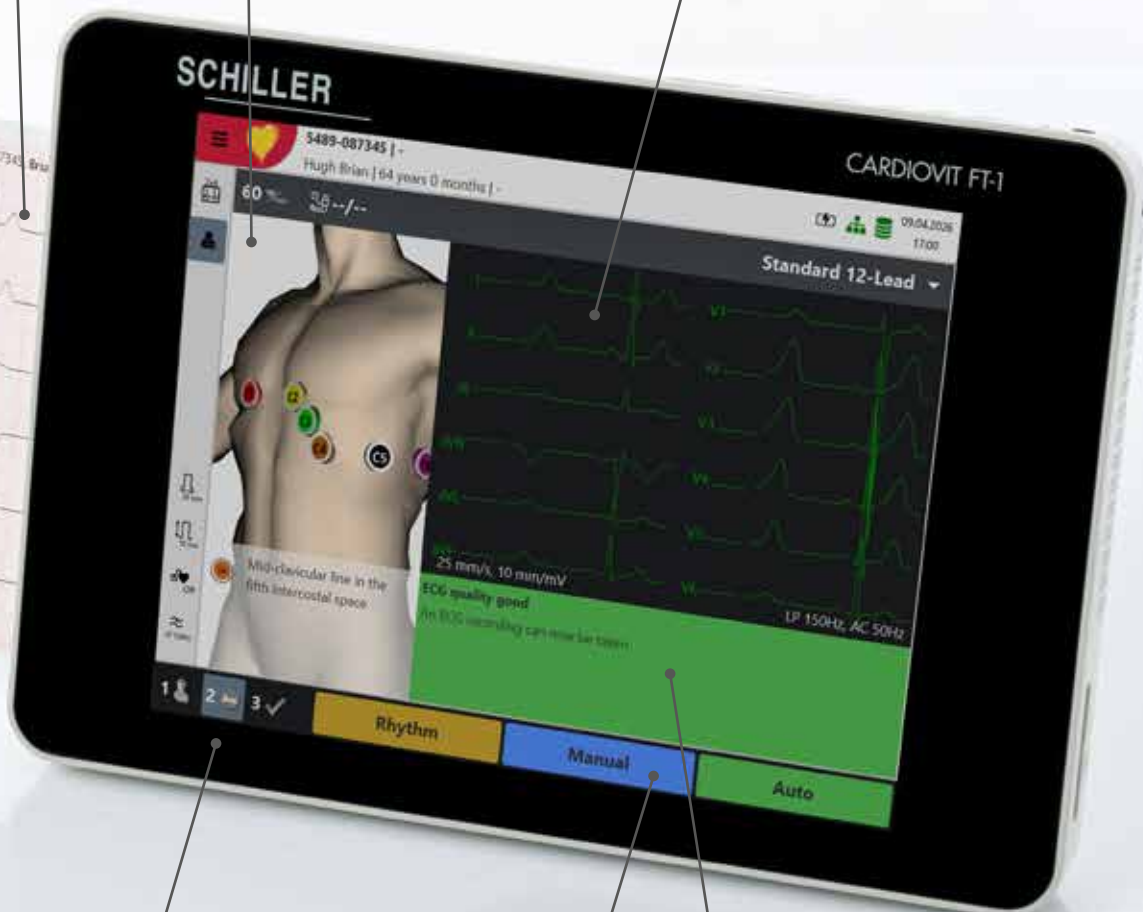

EASY 1-2-3 STEPS

1. Enter or download patient data
 2. Check signal quality and acquire an ECG
 3. Review ECG before storing, exporting or printing
- Users are free to navigate between the steps.

LCD TOUCH SCREEN

8-inch high-resolution multi-touch screen

LEAD QUALITY CHECK

ECG overview of all 12 leads including a colour-coded lead quality check



Americas
SCHILLER Americas Inc.
 Doral, Florida 33172

North America:
 Phone +1 786 845 06 20
 Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com

Latin America & Caribbean
 Phone +1 305 591 11 21
 Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com



Asia
SCHILLER Asia-Pacific / Malaysia
 52200 Kuala Lumpur, Malaysia
 Phone +603 6272 3033
sales@schiller.com.my
www.schiller-asia.com



Austria
SCHILLER Handelsgesellschaft m.b.H.
 A-4040 Linz
 Phone +43 732 70 99 0
sales@schiller.at
www.schiller.at



China
Alfred Schiller (Beijing) Medical Technology Co., Ltd.
 100102 Beijing, China
 Phone +86 10 84565453
info@schillerchina.com
www.schillermedical.cn



Croatia
SCHILLER d.o.o.
 10000 Zagreb
 Phone +385 1 4059 559
info@schiller.hr
www.schiller.hr



France
SCHILLER Médical
 F-67160 Wissembourg
 Phone +33 3 88 63 36 00
info@schiller.fr
www.schiller-medical.fr



France (distribution France)
SCHILLER France S.A.S.
 F-77608 Bussy St Georges
 Phone +33 1 64 66 50 00
contact@schillerfrance.fr
www.schiller-france.com



Germany
SCHILLER Medizintechnik GmbH
 D-85622 Feldkirchen b. München
 Phone +49 89 62 99 81 0
info@schillermed.de
www.schillermed.de



India
SCHILLER Healthcare India Pvt. Ltd.
 Mumbai – 400 059, India
 Phone +91 22 6152 3333 / 2920 9141
sales@schillerindia.com
www.schillerindia.com



Poland
SCHILLER Poland Sp. z o.o.
 PL-02-729 Warszawa
 Phone +48 22 843 20 89 / +48 22 647 35 90
schiller@schiller.pl
www.schiller.pl



Russia & C.I.S.
AO SCHILLER.RU
 119049 Moscow, Russia
 Phone +7 (495) 970 11 33
mail@schiller.ru
www.schiller.ru



Serbia
SCHILLER d.o.o.
 11010 Beograd
 Phone +381 11 39 79 508
info@schiller.rs
www.schiller.rs



Slovenia
SCHILLER d.o.o.
 2310 Slovenska Bistrica
 Phone +386 2 843 00 56
info@schiller.si
www.schiller.si



Spain
SCHILLER ESPAÑA, S.A.
 ES-28232 Las Rozas/Madrid
 Phone +34 91 713 01 76
schiller@schiller.es
www.schiller.es



Switzerland
SCHILLER Schweiz AG
 CH-8912 Obfelden
 Phone +41 44 744 30 00
sales@schiller-schweiz.ch
www.schiller-schweiz.ch



Turkey
SCHILLER TÜRKİYE
 Okmeydanı-Sisli – Istanbul
 Phone +90 212 210 8681
info@schiller.com.tr
www.schiller-turkiye.com



United Kingdom
SCHILLER UK Ltd.
 Bromborough, Wirral,
 Merseyside, CH62 3RJ
 Phone +44 1333 312150
sales@schilleruk.com
service@schilleruk.com
www.schilleruk.com

SCHILLER reserves the right to make changes in design, specifications and features shown or discontinue the product described at any time without notice or obligation. Images may be representative. Some products may not be regulatory cleared or released for sale in all markets. Contact your SCHILLER representative for the most current information and availability.

