

# ERG 911 L ERG 911 LS

Medical couch/semi-couch safety ergometers



**SCHILLER**  
The Art of Diagnostics

# ERG 911 L AND ERG 911 LS

Medical couch ergometer and  
special ergometer for echocardiography during exercise

SCHILLER safety couch/  
semi-couch ergometer:

- ∴ Infinitely electrically  
variable and adjustable
- ∴ Low frame for easy access
- ∴ Quiet and comfortable
- ∴ Configurable



ERG 911 L: couch/semi-couch safety ergometer,  
infinitely variable from 0° - 45°

## ERG 911 L: RECLINED FOR MORE SAFETY

The couch/semi-couch safety ergometer ERG 911 L was developed for use in the cardiology sector. Special safety aspects require the use of a couch/semi-couch safety ergometer, e. g. patients with suspected problems of the cardiovascular system, during exercise following a myocardial infarction, after bypass surgery, or during special use such as cardiac catheter tests. This is when the SCHILLER couch/semi-couch safety ergometer comes into the picture. The use of a couch/semi-couch safety ergometer is also recommended for the examination of elderly or disabled patients. Thanks to the optimal positioning of the patient on the lying surface, a performance identical to bike ergometry is achieved.

## ERG 911 LS: SPECIAL ERGOMETER FOR ECHOCARDIOGRAPHY

The development of the ERG 911 LS is based on the technology of the ERG 911 L couch/semi-couch safety ergometer. Thanks to an additional pivot axis of the couch, the patient can be positioned at an angle of up to 30 degrees to the left during exercise.

This ensures a particularly favourable position for high-quality ultrasound images during stress echocardiography. The stable shoulder support ensures a comfortable and safe position of the patient while the couch is tilted to the left during echocardiography.

Apart from these special applications, the ERG 911 LS can of course still be used as a normal couch/semi-couch safety ergometer.

## ERG 911 LS: INFINITELY VARIABLE TO THE LEFT (0° - 30°)

- ❶ Opening in the lying surface at heart level for convenient ultrasound examinations
- ❷ Stable handlebar and shoulder rest for safe positioning of the patient when the couch is tilted
- ❸ Hip support for increased safety



Two operational modes:

- ❖ Remote operation with a master device (ECG, PC, ...)
- ❖ Standalone operation with either a user-definable step protocol or manual load control

Bluetooth and RS-232 connection

SpO<sub>2</sub> measurement (optional)

Stable, adjustable handlebar

Display rotatable by 180°

Elegant and easy-to-clean casing

High load precision

Very quiet (even at high speed)

Maintenance-free thanks to belt drive

Accurate blood pressure measurement unit with arm rest (optional)

Head rest and seat electrically adjustable

Opening at heart level for stress echocardiography


Infinitely electrically variable to the left (0° - 30°)

Stable steel tube construction for patients weighing up to 160 kg

Infinitely electrically variable (0° - 45°)

Large foot rest with low frame for easy access

Open construction for ergonomic ultrasound examination by the cardiologist

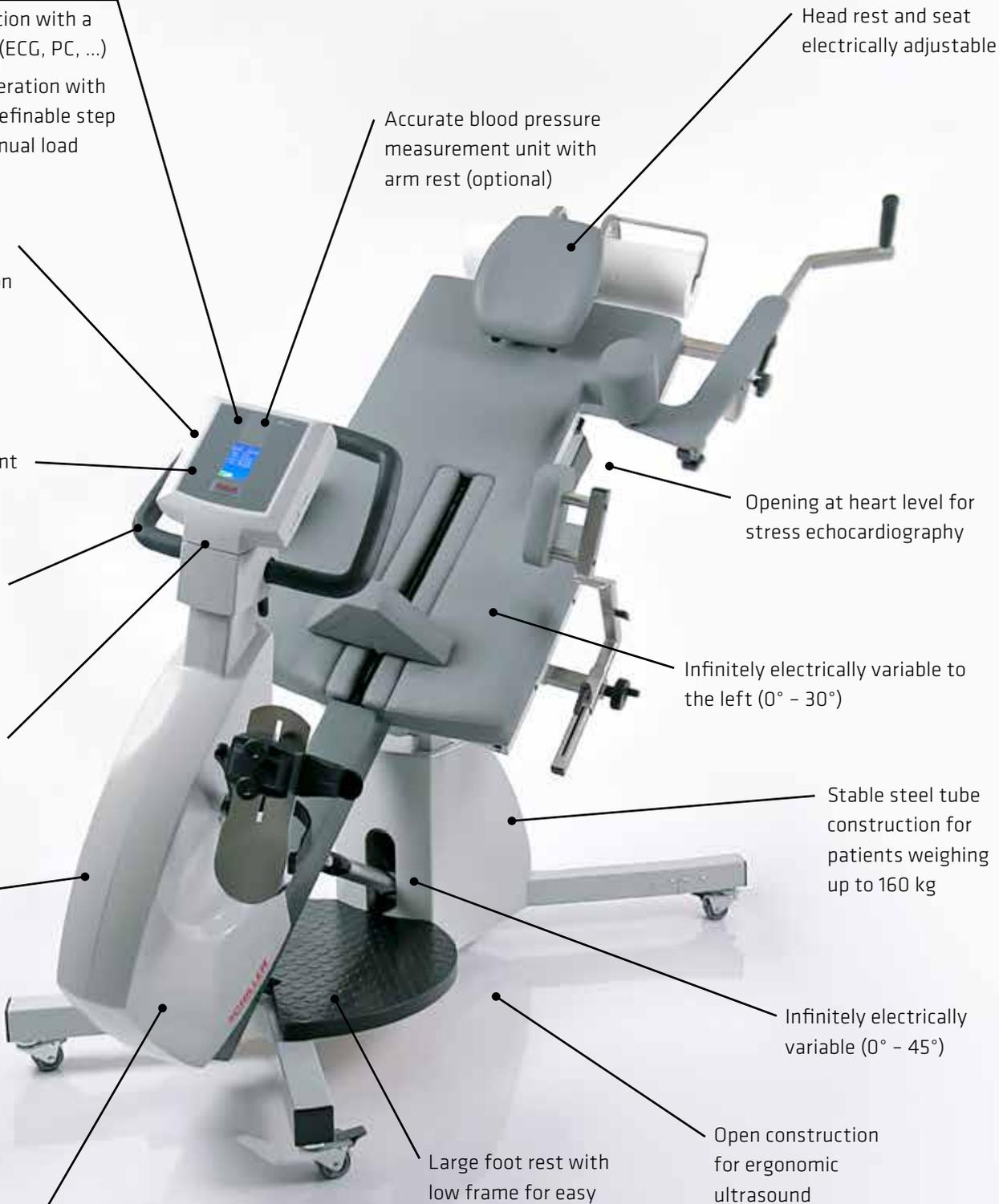



**Americas**  
**SCHILLER Americas Inc.**  
 Doral, Florida 33172

**North America:**  
 Phone +1 786 845 06 20  
 Fax +1 786 845 06 02  
[sales@schilleramericas.com](mailto:sales@schilleramericas.com)  
[www.schilleramericas.com](http://www.schilleramericas.com)

**Latin America & Caribbean**  
 Phone +1 305 591 11 21  
 Fax +1 786 845 06 02  
[sales@schilleramericas.com](mailto:sales@schilleramericas.com)  
[www.schilleramericas.com](http://www.schilleramericas.com)



**Asia**  
**SCHILLER Asia-Pacific / Malaysia**  
 52200 Kuala Lumpur, Malaysia  
 Phone +603 6272 3033  
[sales@schiller.com.my](mailto:sales@schiller.com.my)  
[www.schiller-asia.com](http://www.schiller-asia.com)



**Austria**  
**SCHILLER Handelsgesellschaft m.b.H.**  
 A-4040 Linz  
 Phone +43 732 70 99 0  
[sales@schiller.at](mailto:sales@schiller.at)  
[www.schiller.at](http://www.schiller.at)



**China**  
**Alfred Schiller (Beijing) Medical Technology Co., Ltd.**  
 100102 Beijing, China  
 Phone +86 10 84565453  
[info@schillerchina.com](mailto:info@schillerchina.com)  
[www.schillermedical.cn](http://www.schillermedical.cn)



**Croatia**  
**SCHILLER d.o.o.**  
 10000 Zagreb  
 Phone +385 1 309 66 59  
[info@schillerzg.hr](mailto:info@schillerzg.hr)  
[www.schiller.hr](http://www.schiller.hr)



**France**  
**SCHILLER Médical**  
 F-67160 Wissembourg  
 Phone +33 3 88 63 36 00  
[info@schiller.fr](mailto:info@schiller.fr)  
[www.schiller.fr](http://www.schiller.fr)



**France (distribution France)**  
**SCHILLER France S.A.S.**  
 F-77608 Bussy St Georges  
 Phone +33 1 64 66 50 00  
[contact@schillerfrance.fr](mailto:contact@schillerfrance.fr)  
[www.schiller-france.com](http://www.schiller-france.com)



**Germany**  
**SCHILLER Medizintechnik GmbH**  
 D-85622 Feldkirchen b. München  
 Phone +49 89 62 99 81 0  
[info@schillermed.de](mailto:info@schillermed.de)  
[www.schillermed.de](http://www.schillermed.de)



**India**  
**SCHILLER Healthcare India Pvt. Ltd.**  
 Mumbai - 400 059, India  
 Phone +91 22 6152 3333 / 2920 9141  
[sales@schillerindia.com](mailto:sales@schillerindia.com)  
[www.schillerindia.com](http://www.schillerindia.com)



**Poland**  
**SCHILLER Poland Sp. z o.o.**  
 PL-02-729 Warszawa  
 Phone +48 22 647 35 90  
[schiller@schiller.pl](mailto:schiller@schiller.pl)  
[www.schiller.pl](http://www.schiller.pl)



**Russia & C.I.S.**  
**AO SCHILLER.RU**  
 119049 Moscow, Russia  
 Phone +7 (495) 970 11 33  
[mail@schiller.ru](mailto:mail@schiller.ru)  
[www.schiller.ru](http://www.schiller.ru)



**Serbia**  
**SCHILLER d.o.o.**  
 11010 Beograd  
 Phone +381 11 39 79 508  
[info@schiller.rs](mailto:info@schiller.rs)  
[www.schiller.rs](http://www.schiller.rs)



**Slovenia**  
**SCHILLER d.o.o.**  
 2310 Slovenska Bistrica  
 Phone +386 2 843 00 56  
[info@schiller.si](mailto:info@schiller.si)  
[www.schiller.si](http://www.schiller.si)



**Spain**  
**SCHILLER ESPAÑA, S.A.**  
 ES-28232 Las Rozas/Madrid  
 Phone +34 91 713 01 76  
[schiller@schiller.es](mailto:schiller@schiller.es)  
[www.schiller.es](http://www.schiller.es)



**Switzerland**  
**SCHILLER-Reomed AG**  
 CH-8912 Obfelden  
 Phone +41 44 744 30 00  
[sales@schiller-reomed.ch](mailto:sales@schiller-reomed.ch)  
[www.schiller-reomed.ch](http://www.schiller-reomed.ch)



**Türkiye**  
**SCHILLER TÜRKİYE**  
 Okmeydanı-Sisli - Istanbul  
 Phone +90 212 210 8681 (pbx)  
 Fax +90 212 210 8684  
[info@schiller.com.tr](mailto:info@schiller.com.tr)  
[www.schiller-turkiye.com](http://www.schiller-turkiye.com)



**United Kingdom**  
**SCHILLER UK Ltd.**  
 Anstruther, Fife KY10 3H  
 Phone +44 133 331 2150  
[sales@schilleruk.com](mailto:sales@schilleruk.com)  
[www.schilleruk.com](http://www.schilleruk.com)

SCHILLER reserves the right to make changes in design, specifications and features shown, or discontinue the product described at any time without notice or obligation. Images may be representative. Some products may not be regulatory cleared or released for sale in all markets. Contact your SCHILLER representative for the most current information and availability.

