

medilogFD

The 12-lead Holter ECG for a spatial heart health analysis




SCHILLER
The Art of Diagnostics

medilogFD

The Holter recorder that complements your medilog portfolio:

- ∴ 12-lead Holter ECG allows assessing the full spectrum of cardiac ECG conditions
- ∴ Instantaneous atrial fibrillation/atrial flutter detection based on true P-wave analysis
- ∴ Optional SpO₂ sensor for additional respiratory information

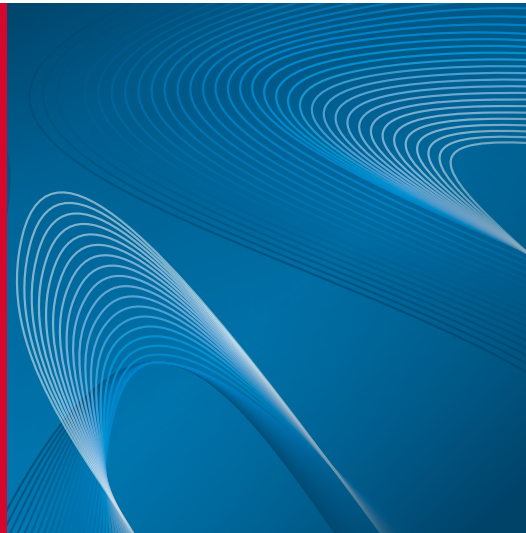

12-LEAD HOLTER ECG FOR A SPATIAL ANALYSIS

The 12-lead Holter ECG is preferred for detailed and spatial heart analysis. Physicians, cardiologists, and electrophysiologists find great value in industry-leading technology:

- ∴ Regional assessment of myocardial ischemia. Distinguish anterior, septal, lateral, inferior region in the dynamic setting of everyday life.
- ∴ Identify the origin of ectopic beats, whether they are from the left or right ventricle.

Superior in identifying abnormal heartbeats that start in the right or left atrium.

MADE TO SUPPRESS ARTIFACTS

The medilogFD Holter recorder features a superior resting ECG-grade 128,000Hz sampling rate on 12 leads. This is used for real-time P-wave analysis, sophisticated artifact suppression, and motion detection, enabling a speedy evaluation.

HOLISTIC APPROACH

Instantaneous detection of atrial fibrillation onset, analysis of respiratory episodes, and high-resolution HRV are just some of the exceptional features of the medilogFD.

DUAL-BATTERY CONCEPT

The high-quality internal battery is recharged via USB or wirelessly on an optional charging station. Moreover, easily available standard AAA batteries can be used. This allows patients to be screened for up to five days on a single charge.

COMPLETES THE MEDILOG PORTFOLIO

medilogFD completes your medilog line-up with medilogAR (3-lead Holter ECG, up to 14 days runtime) and medilog DARWIN2 (high-end ECG analysis software).

High-contrast colour display for easy operation. Live ECG preview: 3 channels shown at a time, all 12 can be scrolled through

Multifactor signal check based on:

- ❖ Impedance
- ❖ Amplitude
- ❖ Amount of noise

Microphone to record patient ID. In an anonymous start, a voice recording is a must to prevent patient mix-up.

Shock-proof thanks to its robust design

Large and easy-to-operate buttons

Visible and tactile position markers for easy patient cable connection

Eco-friendly: runs on an internal rechargeable battery, or, if needed, on normal AAA batteries for fast turnaround times

Data exchange via USB with robust state-of-the-art USB-C connector

Removeable patient cable with 10 wires, tested for more than 4000 connection cycles



Wireless SpO₂ sensor, attached for the night

SCREEN FOR RESPIRATION EPISODES

Thanks to ECG-derived respiration recording (EDR), the medilogFD can screen for potential respiratory episodes during sleep. The optional SpO₂ sensor, connected via Bluetooth, allows for additional respiratory information.



Americas
SCHILLER Americas Inc.
Doral, Florida 33172

North America:
Phone +1 786 845 06 20
Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com

Latin America & Caribbean
Phone +1 305 591 11 21
Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com



France
SCHILLER Médical
F-67160 Wissembourg
Phone +33 3 88 63 36 00
info@schiller.fr
www.schiller.fr



France (distribution France)
SCHILLER France S.A.S.
F-77608 Bussy St Georges
Phone +33 1 64 66 50 00
contact@schillerfrance.fr
www.schiller-france.com



Serbia
SCHILLER d.o.o.
11010 Beograd
Phone +381 11 39 79 508
info@schiller.rs
www.schiller.rs



Slovenia
SCHILLER d.o.o.
2310 Slovenska Bistrica
Phone +386 2 843 00 56
info@schiller.si
www.schiller.si



Asia
SCHILLER Asia-Pacific / Malaysia
52200 Kuala Lumpur, Malaysia
Phone +603 6272 3033
sales@schiller.com.my
www.schiller-asia.com



Germany
SCHILLER Medizintechnik GmbH
D-85622 Feldkirchen b. München
Phone +49 89 62 99 81 0
info@schillermed.de
www.schillermed.de



Spain
SCHILLER ESPAÑA, S.A.
ES-28232 Las Rozas/Madrid
Phone +34 91 713 01 76
schiller@schiller.es
www.schiller.es



Austria
SCHILLER Handelsgesellschaft m.b.H.
A-4040 Linz
Phone +43 732 70 99 0
sales@schiller.at
www.schiller.at



India
SCHILLER Healthcare India Pvt. Ltd.
Mumbai – 400 059, India
Phone +91 22 6152 3333 / 2920 9141
sales@schillerindia.com
www.schillerindia.com



Switzerland
SCHILLER Schweiz AG
CH-8912 Obfelden
Phone +41 44 744 30 00
sales@schiller-schweiz.ch
www.schiller-schweiz.ch



China
Alfred Schiller (Beijing) Medical Technology
Co., Ltd.
100102 Beijing, China
Phone +86 10 84565453
info@schillerchina.com
www.schillermedical.cn



Poland
SCHILLER Poland Sp. z o.o.
PL-02-729 Warszawa
Phone +48 22 843 20 89 / +48 22 647 35 90
schiller@schiller.pl
www.schiller.pl



Turkey
SCHILLER TÜRKIYE
Okmeydani-Sisli – Istanbul
Phone +90 212 210 8681 (pbx)
Fax +90 212 210 8684
info@schiller.com.tr
www.schiller-turkiye.com



Croatia
SCHILLER d.o.o.
10000 Zagreb
Phone +385 1 309 66 59
info@schiller.hr
www.schiller.hr



Russia & C.I.S.
AO SCHILLER.RU
119049 Moscow, Russia
Phone +7 (495) 970 11 33
mail@schiller.ru
www.schiller.ru



United Kingdom
SCHILLER UK Ltd.
Anstruther, Fife KY10 3H
Phone +44 1333 312150
sales@schilleruk.com
www.schilleruk.com

SCHILLER reserves the right to make changes in design, specifications and features shown, or discontinue the product described at any time without notice or obligation. Images may be representative. The illustration of patient cables and screen contents show the system using an IEC cable. Depending on your region, the unit is supplied with an AHA patient cable. Some products may not be regulatory cleared or released for sale in all markets. Contact your SCHILLER representative for the most current information and availability.

