

EASY PULSE®

Efficient and compact mechanical CPR
thanks to technological innovation


Latest Generation of
Automatic CPR



SCHILLER
The Art of Diagnostics

EASY PULSE®

SCHILLER has taken the next step in improving CPR quality further with the launch of the unbelievably small and light mechanical chest compression device: the EASY PULSE.

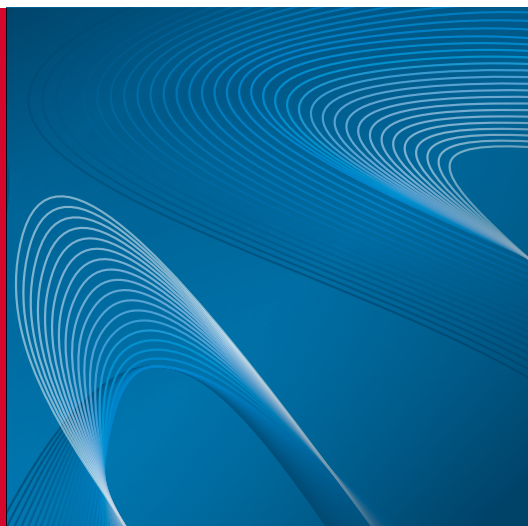

Performing manual chest compressions well for an extended period of time is almost impossible. Not only is it physically demanding, but other actions are also required simultaneously, such as vital signs monitoring.

SCHILLER's EASY PULSE is the solution for more effective resuscitation: this portable, stand-alone device delivers chest compressions automatically at a consistent rate and depth. It is directly attached to the patient's upper body and can thus be used in any situation, regardless of ambient conditions. Especially during transport, a high-quality uninterrupted chest compression can be guaranteed.

The EASY PULSE is such a simple device that a single person can operate it. Its size and weight are unparalleled: the EASY PULSE weighs 3.5 kg, opening up entirely new possibilities e.g. for air rescue services.

LATEST GENERATION OF AUTOMATIC CPR DEVICES

The EASY PULSE has a unique combination of stamp and band which allows a 3D compression and therefore maximum efficiency.

The weight of the device is evenly spread over the baseplate and not on a single point and therefore has no influence on chest compressions nor on the recoil.



SLIDER

The slider enables a relieved preparation of the patient for mechanical chest compressions. It also takes very little time to attach the device thanks to the practical buckle system – less than 15 seconds for trained personnel.



CHARGER

The rechargeable battery can be charged in a dual charger (less than 2 hours) or in turbo mode (less than 50 minutes) and lasts for approx. 45 minutes.



EASE OF USE

In an emergency every move has to be right. The EASY PULSE is simple and intuitive to operate: pressing two buttons is enough to initiate resuscitation immediately.



EXTERNAL POWER SUPPLY AND DIRECT CONNECTION CABLE

A power supply unit for 100VAC to 240VAC input voltage is supplied. The EASY PULSE can be directly operated using a 12VDC - 24VDC power source via a special connection cable.

TECHNOLOGICAL INNOVATIONS

Lightweight

- ❖ Unbelievably small and light electrically driven CPR device, weighing 3.5 kg
- ❖ Ideal for limited spaces

Comprehensive

- ❖ Multidirectional chest compressions (3D)
- ❖ High-quality chest compressions of consistent rate and depth, which is impossible to achieve with manual compressions
- ❖ Possibility of 30:2 compression-ventilation cycle

Simple

- ❖ Thanks to slider and buckle system, it can be easily attached to the upper body
- ❖ Easy to position – easy to operate

Additional features

- ❖ Replaceable battery with charge level indicator
- ❖ Autonomous operation for 45 minutes (with one battery)
- ❖ Connections: external DC input, USB



Americas
SCHILLER Americas Inc.
Doral, Florida 33172

North America:
Phone +1 786 845 06 20
Fax +1 786 845 06 02
sales@schilleramerica.com
www.schilleramerica.com

Latin America & Caribbean
Phone +1 305 591 11 21
Fax +1 786 845 06 02
sales@schillerla.ch
www.schillerla.ch



Asia
SCHILLER Asia-Pacific / Malaysia
52200 Kuala Lumpur, Malaysia
Phone +603 6272 3033
Fax +603 6272 2030
sales@schiller.com.my
www.schiller-asia.com



Austria
SCHILLER Handelsgesellschaft m.b.H.
A-4040 Linz
Phone +43 732 70 99
Fax +43 732 757 000
sales@schiller.at
www.schiller.at



China
Alfred Schiller (Beijing) Medical Equipment
Co., Ltd.
100015 Beijing, China
Phone +86 010 52007020
Fax +86 010 52007020-8016
info@schillerchina.com
www.schillermedical.cn



Croatia
Schiller medicinski instrumenti d.o.o.
10000 Zagreb
Phone +385 1 309 66 59
Fax +385 1 309 66 60
info@schillerzg.hr
www.schiller.hr



France
SCHILLER Médical S.A.S.
F-67162 Wissembourg
Phone +33 3 88 63 36 00
Fax +33 3 88 63 36 49
info@schiller.fr
www.schiller-medical.com



France (distribution France)
SCHILLER France S.A.S.
F-77608 Bussy St Georges
Phone +33 1 64 66 50 00
Fax +33 1 64 66 50 10
contact@schillerfrance.fr
www.schiller-france.com



Germany
SCHILLER Medizintechnik GmbH
D-85622 Feldkirchen b. München
Phone +49 89 62 99 81 0
Fax +49 89 609 50 90
info@schillermed.de
www.schillermed.de



Hungary
SCHILLER Diamed Ltd.
H-1141 Budapest
Phone +36 1 383 4780 / +36 1 460 9491
Fax +36 1 383 4778
info@schillerhungary.hu
www.schillerhungary.hu



India
SCHILLER Healthcare India Pvt. Ltd.
Mumbai - 400 059, India
Phone +91 22 6152 3333 / 2920 9141
Fax +91 22 2920 9142
sales@schillerindia.com
www.schillerindia.com



Poland
SCHILLER Poland Sp. z o.o.
PL-02-729 Warszawa
Phone +48 22 843 20 89 / +48 22 647 35 90
Fax +48 22 843 20 89
schiller@schiller.pl
www.schiller.pl



Russia & C.I.S.
AO SCHILLER.RU
125124 Moscow, Russia
Phone +7 (495) 970 11 33
Fax +7 (495) 956 29 10
mail@schiller.ru
www.schiller.ru
www.schiller-cis.com



Serbia
SCHILLER d.o.o.
11010 Beograd
Phone +381 11 39 79 508
Fax +381 11 39 79 518
info@schiller.rs
www.schiller.rs



Slovenia
SCHILLER d.o.o.
2310 Slovenska Bistrica
Phone +386 2 843 00 56
Fax +386 2 843 00 57
info@schiller.si
www.schiller.si



Spain
SCHILLER ESPAÑA, S.A.
ES-28232 Las Rozas/Madrid
Phone +34 91 713 01 76
Fax +34 91 355 79 33
schiller@schiller.es
www.schiller.es



Switzerland
SCHILLER-Reomed AG
CH-8953 Dietikon
Phone +41 44 744 30 00
Fax +41 44 740 37 10
sales@schiller-reomed.ch
www.schiller-reomed.ch



Turkey
SCHILLER TÜRKİYE
Okmeydanı-Sisli - Istanbul
Phone +90 212 210 8681 (pbx)
Fax +90 212 210 8684
info@schiller.com.tr
www.schiller-turkiye.com



United Kingdom
SCHILLER UK
Bellshill, ML4 3PR
Phone +44 1698 744 505
Fax +44 1698 744 474
sales@schilleruk.com
www.schilleruk.com

