


CARDIOVIT AT-102 G2

Easy, smart and state-of-the-art

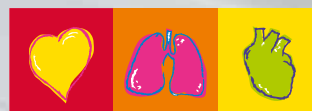

SCHILLER
The Art of Diagnostics

CARDIOVIT AT-102 G2

The new generation ECG
CARDIOVIT AT-102 G2:

- ∴ Robust and reliable for hospitals and clinics
- ∴ Fast bidirectional Wi-Fi communication with HIS/EMR
- ∴ Best signal quality and diagnostic algorithms for adult and paediatric patients



STATE-OF-THE-ART TECHNOLOGY AND CONNECTIVITY

The CARDIOVIT AT-102 G2 features an 8" colour display with touch function keys. Moreover, its sealed alphanumeric keyboard ensures easy cleaning, making it highly suitable for the daily clinical routine. Smart Battery technology guarantees more than 8 hours of ECG recording.

Bidirectional communication allows for easy data access and fast transmission of ECG reports to EMR/HIS systems. Standard Wi-Fi and LAN connectivity, combined with ECG review on display and a large memory, supports paperless workflows and cost saving. When a paper report is necessary, customizable printout are available on the internal high resolution A4 thermal printer.

HIGH EFFICIENCY IN ANY WORKFLOW

ECG acquisition has never been easier and faster: Enter or download patient data with PDQ or worklist, acquire ECG, save, print and transmit data in twenty seconds using only one navigation key.



EASY 1-2-3 STEPS

- 1: Enter or download patient data
- 2: Control signal quality and lead reversal, and perform ECG
- 3: Review ECG before storing, exporting or printing

BEST SIGNAL QUALITY AND DIAGNOSTIC TOOLS

SCHILLER's most advanced algorithms are implemented in the CARDIOVIT AT-102 G2:

PAEDIATRIC ECG: Thanks to very high sampling frequency and a large frequency band up to 250 Hz, the CARDIOVIT AT-102 G2 offers optimal signal quality for paediatric ECG.

RESTING RHYTHM: With the CARDIOVIT AT-102 G2, resting rhythm can be recorded for up to 10 minutes. Full disclosure, HR trend and events can be reviewed, stored and printed, as well as transmitted to HIS/EMR in PDF, DICOM or HL7 format.

Internal high resolution
A4 thermal printer

Sealed alphanumeric keyboard,
easy to clean and highly hygienic

Touch function keys

Navigation keys allow high
efficiency in any workflow

Organize and safeguard
hospital assets and personal
items with a bin

Manoeuvre smoothly and
comfortably with easy-glide
wheels and a handle positioned
at just the right height



EXERCISE ECG: The CARDIOVIT AT-102 G2 can optionally perform basic exercise ECG with a large selection of protocols for treadmills and ergometers. A complete review and report of the test with trends, step table, interpretation template and event list are available. Exercise data can also be transmitted to HIS/EMR in different formats.

ETM WITH ETM SPORT: One of the industry's most renowned, simultaneous 12-lead ECG interpretation program for uncompromising quality and reliability.

ETM Sport is the first automated ECG interpretation module designed for athletes, based on the Seattle Criteria.

CCAA: SCHILLER's Culprit Coronary Artery Algorithm is used in emergency departments and outside the hospital whenever a fast decision on chest pain is required immediately. It locates the obstruction site in the coronary artery and thereby the size of the cardiac area at risk.



Americas
SCHILLER Americas Inc.
 Doral, Florida 33172

North America:
 Phone +1 786 845 06 20
 Fax +1 786 845 06 02
sales@schilleramerica.com
www.schilleramerica.com

Latin America & Caribbean
 Phone +1 305 591 11 21
 Fax +1 786 845 06 02
sales@schillerla.ch
www.schillerla.ch



Asia
SCHILLER Asia-Pacific / Malaysia
 52200 Kuala Lumpur, Malaysia
 Phone +603 6272 3033
 Fax +603 6272 2030
sales@schiller.com.my
www.schiller-asia.com



Austria
SCHILLER Handelsgesellschaft m.b.H.
 A-4040 Linz
 Phone +43 732 70 99
 Fax +43 732 757 000
sales@schiller.at
www.schiller.at



China
Alfred Schiller (Beijing) Medical Equipment Co.,Ltd.
 100015 Beijing, China
 Phone +86 010 52007020
 Fax +86 010 52007020-8016
info@schillerchina.com
www.schillermedical.cn



Croatia
SCHILLER d.o.o.
 10000 Zagreb
 Phone +385 1 309 66 59
 Fax +385 1 309 66 60
info@schillerzg.hr
www.schiller.hr



France
SCHILLER Médical
 F-67160 Wissembourg
 Phone +33 3 88 63 36 00
 Fax +33 3 88 63 36 49
info@schiller.fr
www.schiller-medical.com



France (distribution France)
SCHILLER France S.A.S.
 F-77608 Bussy St Georges
 Phone +33 1 64 66 50 00
 Fax +33 1 64 66 50 10
contact@schillerfrance.fr
www.schiller-france.com



Germany
SCHILLER Medizintechnik GmbH
 D-85622 Feldkirchen b. München
 Phone +49 89 62 99 81 0
 Fax +49 89 609 50 90
info@schillermed.de
www.schillermed.de



Hungary
SCHILLER Diamed Ltd.
 H-1141 Budapest
 Phone +36 1 383 4780 / +36 1 460 9491
 Fax +36 1 383 4778
info@schillerhungary.hu
www.schillerhungary.hu



India
SCHILLER Healthcare India Pvt. Ltd.
 Mumbai - 400 059, India
 Phone +91 22 6152 3333 / 2920 9141
 Fax +91 22 2920 9142
sales@schillerindia.com
www.schillerindia.com



Poland
SCHILLER Poland Sp. z o.o.
 PL-02-729 Warszawa
 Phone +48 22 843 20 89 / +48 22 647 35 90
 Fax +48 22 843 20 89
schiller@schiller.pl
www.schiller.pl



Russia & C.I.S.
AO SCHILLER.RU
 125040 Moscow, Russia
 Phone +7 (495) 970 11 33
 Fax +7 (495) 956 29 10
mail@schiller.ru
www.schiller.ru



Serbia
SCHILLER d.o.o.
 11010 Beograd
 Phone +381 11 39 79 508
 Fax +381 11 39 79 518
info@schiller.rs
www.schiller.rs



Slovenia
SCHILLER d.o.o.
 2310 Slovenska Bistrica
 Phone +386 2 843 00 56
 Fax +386 2 843 00 57
info@schiller.si
www.schiller.si



Spain
SCHILLER ESPAÑA, S.A.
 ES-28232 Las Rozas/Madrid
 Phone +34 91 713 01 76
 Fax +34 91 355 79 33
schiller@schiller.es
www.schiller.es



Switzerland
SCHILLER-Reomed AG
 CH-8953 Dietikon
 Phone +41 44 744 30 00
 Fax +41 44 740 37 10
sales@schiller-reomed.ch
www.schiller-reomed.ch



Turkey
SCHILLER TÜRKİYE
 Okmeydanı-Sisli - Istanbul
 Phone +90 212 210 8681 (pbx)
 Fax +90 212 210 8684
info@schiller.com.tr
www.schiller-turkiye.com



United Kingdom
SCHILLER UK
 Bellshill, ML4 3PR
 Phone +44 1698 744 505
 Fax +44 1698 744 474
sales@schilleruk.com
www.schilleruk.com

Article no.: 2.500767
 2018/04, EN



SCHILLER
 The Art of Diagnostics